

SUCTION MASTER  
PURISUC



유해한 액체나 폐시료를 신속히  
제거하기 위한 컴팩트한 석션 펌프

STANDARD SPECIFICATION

Model	PURISUC 2500ml	
Size	Bottle(W)	150Ø
	Outside(W x D x H)mm	200 x 340 x 226
	Weight(Kg)	18
Controller	Manual On/Off switch, Foot switch	
Outside Material	ABS	
Bottom Material	Epoxy powder coated steel	
Remote Operation	Foot switch connected electricity cable	
Suction Pump	Oilless pump, 600mmHg, 6Kgf/cm , 12L/min	
Power	230V ±10%, 50-60Hz, 1Phase	
Power Consumption	105W	
Safety	Fuse installed noise filter, Water level sensor in bottle	

FEATURE PLUS

컴팩트한 사이즈  
클린벤치 하부에서의 사용에 적합

Application Diagram



## MAIN FEATURE




### 수위센서 장착으로 사용안전성 증대

수위센서가 장착되어 있어서 용기 내 내용물이 일정 높이에 도달하면, 자동으로 흡입기능이 정지하도록 디자인 되어 있습니다.



### 핸즈프리 풋 스위치 기능으로 사용 편의 증대

풋 패드 타입의 On/off 스위치를 통하여, 사용자가 손을 사용하지 않고 쉽게 작동 가능합니다.

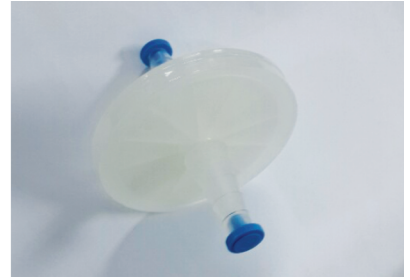

### 역류방지 필터 장착으로 내구성 증대

역류방지 필터 장착으로 내구성 증대 역압에 의한 내용물의 펌프 방향으로의 역류를 차단하여 펌프의 고장 위험을 제거하였기에 전반적인 장비 내구성이 강화되었습니다.

## ORDER INFORMATION

구분	ITEM	ORDER CODE
		PURISUC 2500
제품	Main Body	<b>PUS-25-MB</b>
스페어	Bottle(2500ml) Only	PUS-25-SP1
	Bottle(2500ml)+Cap Set	PUS-25-SP2
	Foot Switch	PUS-25-SP3
	Water Level Sensor	PUS-25-SP4